

Software License Set Up


1. Run `install.bat` as administrator (right click > "Run as administrator")
2. The software license is managed via the Admin Control Center (ACC): a browser-based UI. Access the ACC (on port 1947) in the browser (`localhost:1947`).
3. Once in the ACC, go to the [Sentinel Keys](#) page, and get a C2V file by clicking the *Fingerprint* button in the *Actions* column. Send the C2V file to `support@voronoihealthanalytics.com`.

Sentinel Keys

Location	Vendor	Key ID	Key Type	Configuration	Version	Sessions	Actions
Local	32088 (32088)		Reserved for New SL Key	SL	8.23		Fingerprint

4. After sending the C2V file to Voronoi Health Analytics, you will be sent a V2C file. The V2C file contains the information required to update the software license so that it contains the license for the DAFS application.

Apply the V2C file in the [Update/Attach](#) page of the ACC.

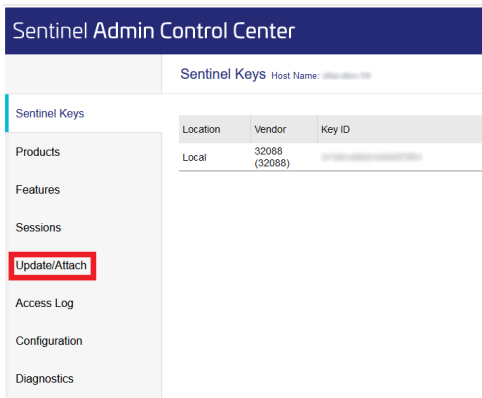

5. In addition to the V2C file, you will receive a `*VHALm.exe` file. Run `*VHALm.exe`.

Additional Information

- For license updates a C2V file can be downloaded by clicking the C2V button on the [Sentinel Keys](#) page of the ACC.

Sentinel Keys

Location	Vendor	Key ID	Key Type	Configuration	Version	Sessions	Actions
Local	32088 (32088)	#####	HASP SL AdminMode		8.23		Products Features Sessions Certificates C2V

- To remove the software license run `uninstall.bat` as administrator (right click > "Run as administrator")